





1.
choose temperatur



2.
switch heat on (takes 3-5 min.)

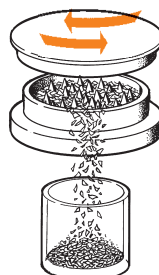


3.
grind substance

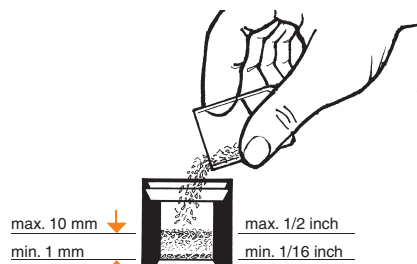

max. 10 mm
min. 1 mm
max. 1/2 inch
min. 1/16 inch

4.
fill filling chamber

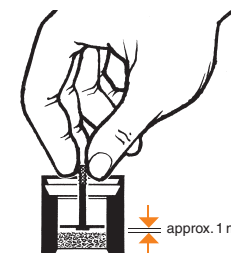

5.
insert filling chamber insert


approx. 3 sec




6.
preheat substance if
control lamp is off




7.
„click“ valve balloon
on filling chamber

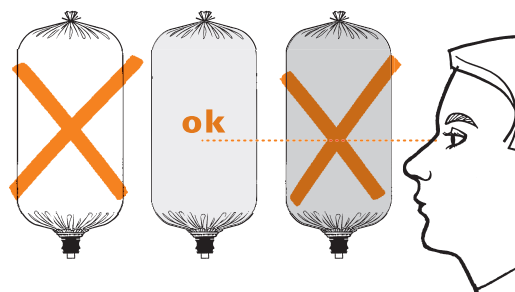

8.
watch vapor density

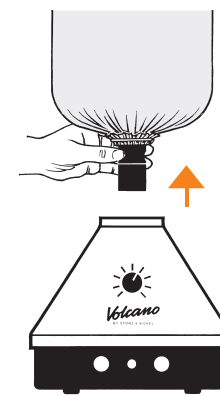


9.
remove full balloon unit

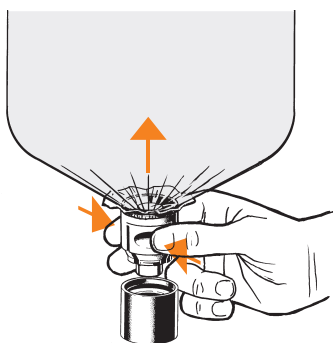


10.
separate valve balloon
from filling chamber




11.
Caution! Causes burns!
don't leave filling chamber on VOLCANO

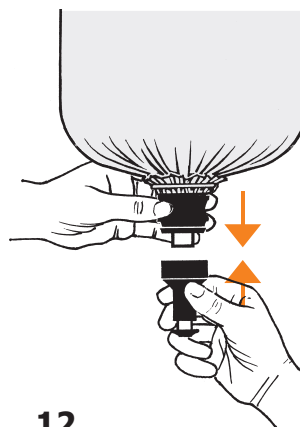

12.
„click“ mouthpiece on
valve balloon

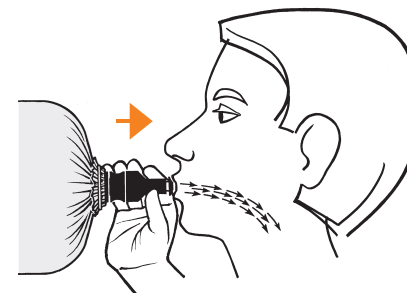


13.
press mouthpiece against
lips to inhale

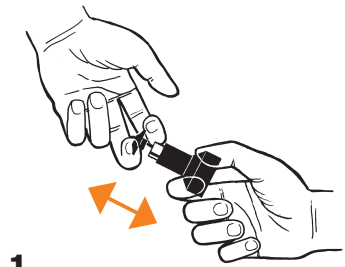


1.
draw of lip piece



2.
dismantled
mouthpiece

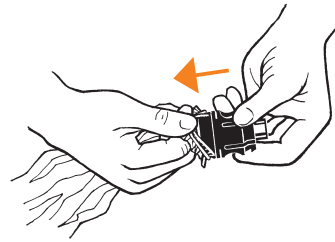

3.
draw of
balloon

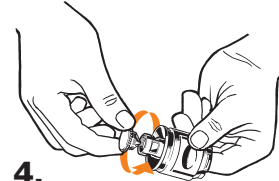

4.
srew off
valve cap




5.
dismantled
valve




6.
screw off nut from
filling chamber insert



7.
dismantled filling
chamber insert

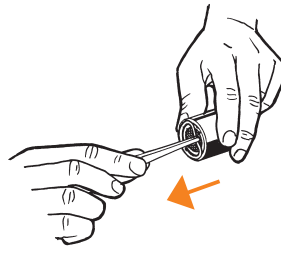

8.
draw off filling
chamber ring

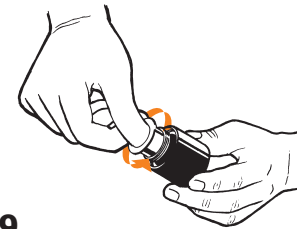


9.
srew off filling
chamber cylinder




10.
dismantled
filling chamber




11.
manual cleaning
with alcohol (min 90 %)

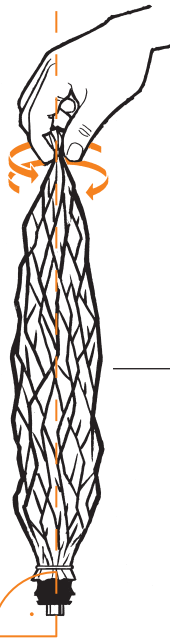

12.
cut off balloon

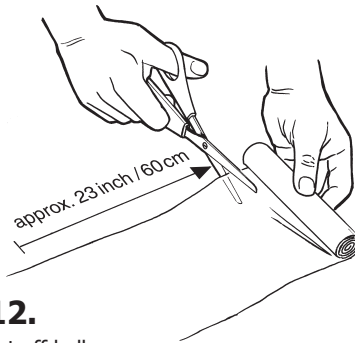

13.
tie up balloon

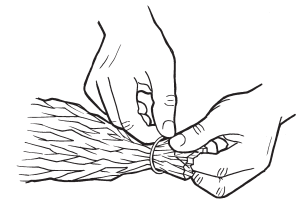


14.
attach first balloon
fixation ring




15.
attach balloon to valve


16.
align balloon to valve

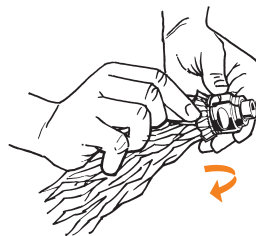

17.
evert balloon mantle

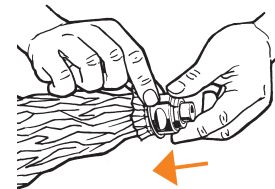

18.
attach second balloon
fixation ring to secure

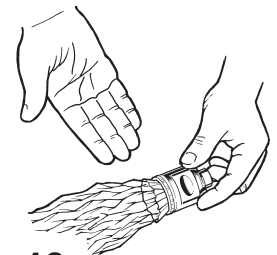

19.
ready to use